
ARDESTO

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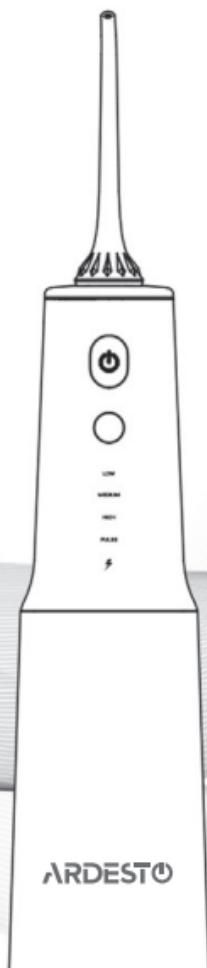HEALTHY TEETH AND GUMS! EASY CLEANING OF ORTHODONTIC CONSTRUCTIONS!


You can effectively and easily take care of your gums, teeth, and beauty of your smile, using Ardesto irrigator on a daily basis! With the combination of water pressure and pulsation, the irrigator removes harmful bacteria and food debris deep between teeth and below gumline braces and other orthodontic constructions where brushing and even dental floss won't reach.

Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria



Flossing has limited access below the gumline and may damage gums

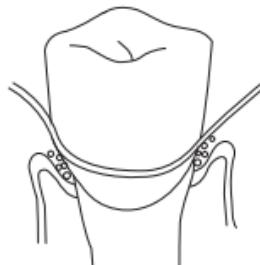Irrigator Clean

Removes trapped food debris and harmful bacteria with high frequency pulse water floss.


If you use irrigator on a daily basis, it gives you fresh breath, makes gums and teeth healthier, and helps to clean orthodontic constructions! Consult with your dentist on the irrigator modes usage.

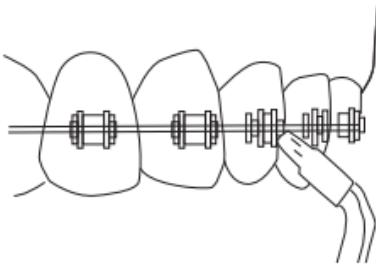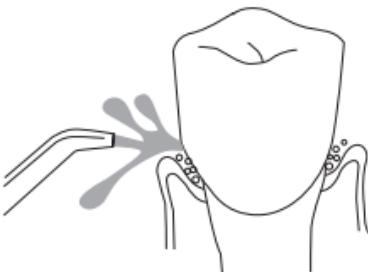CONTENTS

| | |
|--------------------------------------|----|
| Product Description | 15 |
| Operation Steps..... | 16 |
| Using Tips | 16 |
| Disassembly & Assembly Examples..... | 17 |
| Equipment | 20 |
| Product Accessories | 20 |
| UL Standard | 21 |
| Battery Care and Disposal | 22 |
| Common issues | 22 |
| Warranty Conditions | 24 |

PRODUCT DESCRIPTION

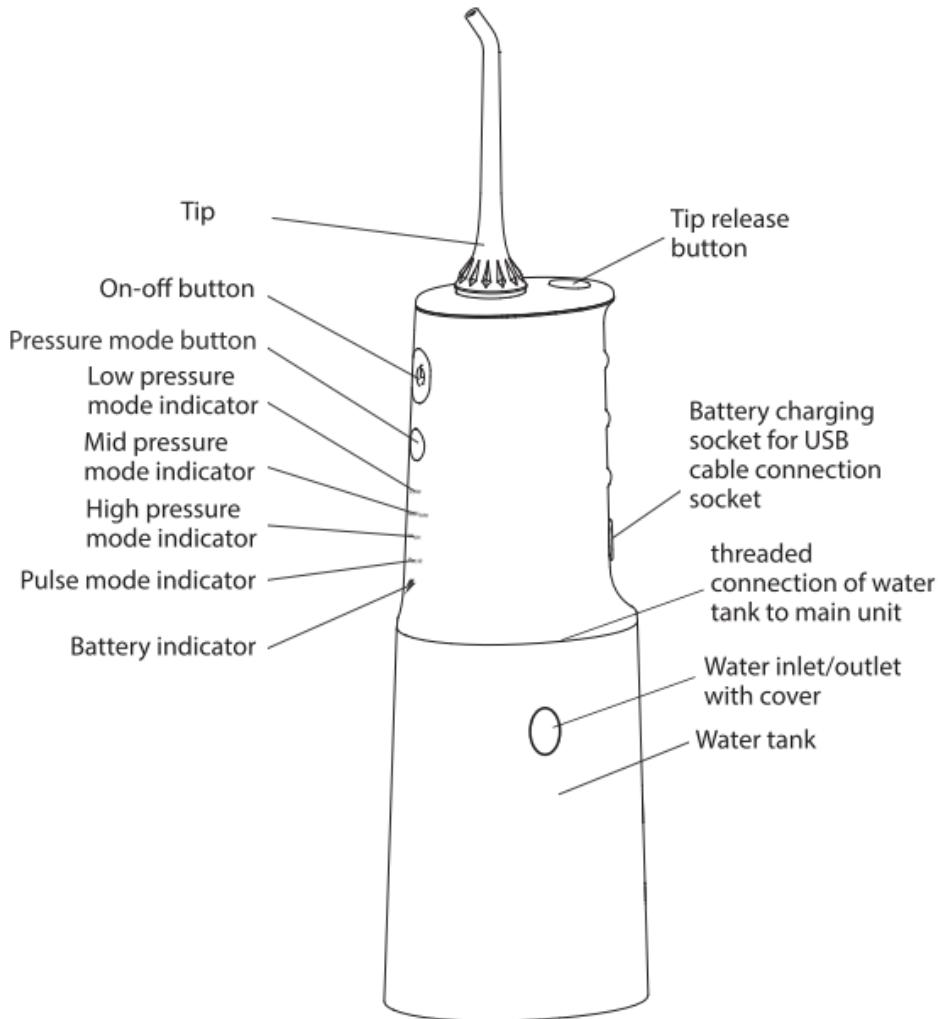
OPERATING STEPS

Start

IMPORTANT! Please fully charge the Water flosser for the first use by plugging the USB cable into the charging socket with power. The product cannot be used while charging.

USB adapter is not included. The USB cable is compatible with any 220V/5V, 1A USB power adapter.

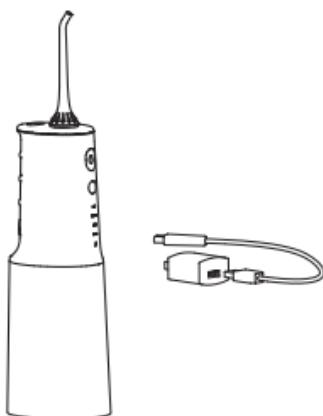
USING TIPS

Tilt down to place the tip in the mouth, and point the tip at the tooth, then turn on the device.

Make water jet at 90° to the gingiva, and close your lips slightly to avoid water spray. For best results, start with the molars (back-teeth) and clean towards the front teeth. Slide the tip along the gingiva line and rest briefly between teeth until it is cleaned around gingiva and between teeth.

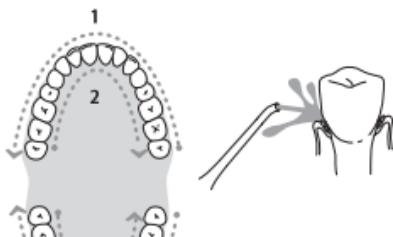
DISASSEMBLY & ASSEMBLY EXAMPLES


After-use

Drainage of residual liquid

Turn off the device, and open the back tank lid with the injection hole facing down (shown in Fig. 1), then turn on the device after clearing in the internal channels of the irrigator the residual liquid in the tank, to drain the residual liquid in the until there is no more water to drain, then turn off the device.

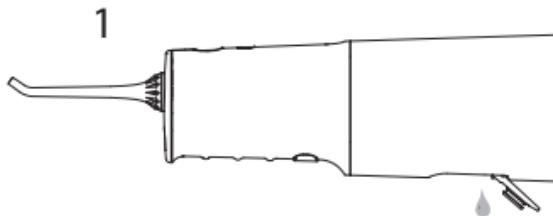Separate water tank from the main unit of the irrigator to wash or dry it. It's easy to do. Hold the product body with your left hand and hold the tank with the right hand, as shown in Fig. 2, and turn it counterclockwise to pull down the water tank, as shown in Fig. 3.

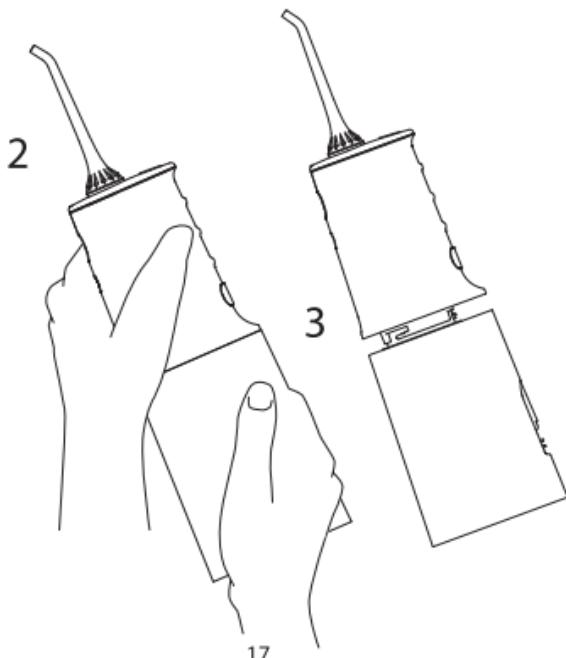
Installation and removal of tips


Insert the nozzle into the middle hole of the tip knob until it clicks.

As shown in Fig. 4, open the tank lid and fill water into the tank;
Note: Make sure that the device is off while filling the tank.

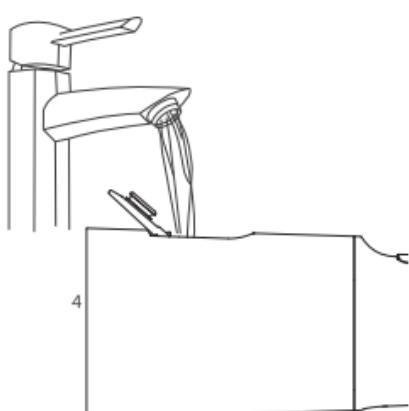If the tip needs replacing or removing, press the release button and while remove the tip.

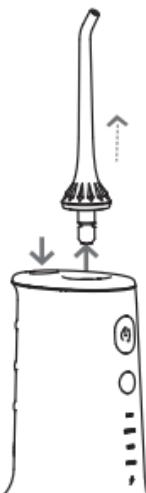
- This product has memory function, so default work as the mode of last using.
- Press the mode button in the power-on state to switch between the low, mid, high and pulse modes; the working gear position and the power display are always on, and the other non-working gears are off.
- Mode selection can also be made when power off, and you can choose your favorite mode and then turn it on. Running in any mode for more than 120 seconds, the product will shut down automatically.
- When the product is charging, the power indicator light flashes; when the power indicator light is on, the product is fully charged.
- During use, when the power indicator light changes from white to red, the product is low battery; when the red light is flashing, the product needs charging.
- This product cannot be used while charging.
- Rotate the tip knob to adjust the direction of the water flow.
- Because gingiva bearing capacity differs from people to people, the first use of the Water Flosser may be accompanied with bleeding gingiva. This is normal. Please feel free to use it!
- In the first use of the portable water flosser, please start from the lowest mode and gradually shift to the higher with oral adaptation.

EQUIPMENT

Main body – 1 pcs
Water tank – 1 pcs
Jet tips – 2 pcs
USB cable – 1 pcs
User manual with warranty card – 1 pcs

PRODUCT ACCESSORIES

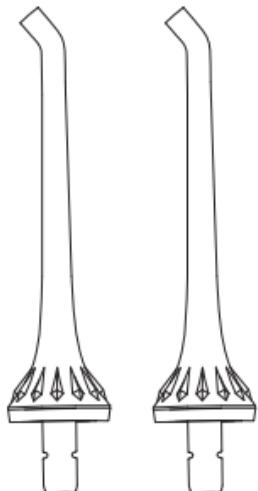Tips x 2

USB cable x 1

Replaceable tip/accessories are available from retail stores

UL STANDARD

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.

DANGERS:

- To reduce the risk of electrocution:
 - Do not handle charger with wet hands.
 - Do not immerse in water or other liquid.
 - Do not use while bathing.
 - Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
 - Do not reach for a product that has fallen into water. Unplug immediately.
 - Check the charger cord for damage before the first use and during the life of the product.

To reduce the risk of burns, electric shock, fire, serious injury or damage issues:

Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.

Do not use this product if it has a broken cord or plug. In case of not working properly, or it has been dropped into any liquid, contact your local distributor.

Do not direct water flow under the tongue, into the ear, nose or other delicate organs. This product is capable of producing pressures that may cause serious damage to these organs. Refer to operating instructions for correct usage.

Use this product only as indicated in these instructions or suggestions by your dental professional.

It's normal phenomenon if gum bleed a little when using water flosser for the first time. If bleed much, please stop using and ask for your dentist advice.

Fill reservoir with water or other solutions recommended by dental professional.

Do not drop or insert any foreign object into the product.
Keep charger and cord away from heated surfaces.
Do not operate where oxygen or aerosol sprays are being used.
Do not use mouth washes/rinses that contain iodine, bleach, or tea treeoil.
Do not use while wearing any oral jewelry. Please put off prior to use.
Do not use if you have an open wound on your tongue or in your mouth.
If your physician advise you to receive antibiotic premedication,
you should consult your dentist before using this instrument or any
other oral hygiene aid.
For household use only. Do not use outdoors.
Closely instruct and supervise children and individualswith special
needs in the proper use of this product.
Children should be supervised to ensure that they donot play with
the appliance.

BATTERY CARE AND DISPOSAL

Battery Maintenance

- Lithium-ion battery: 3.7V, 1900mAh One 4-hour charge will be enough for a week. This type of battery is allowed for air transportation in your checked-in or cabin baggage. Please contact your air carrier for the transportation conditions.
- The device must be fully charged for 24 hours for the first use, since then charging about 4 hours.
- Plug the USB cable into the charging socket and power on;
- If it is idle for a long time, fully charge the device;
- Charge immediately once low-battery warning so that to utilize battery well.

Disposing of the Battery

The device is equipped with a lithium battery, which will be recycled at the end of life.

Please take batteries, which are non-replaceable, to your local recycling center.

FREQUENTLY ASKED QUESTIONS

| Issues | Reason | Treatment |
|---|---|---|
| Leakage between the tip and the station | Broke seal ring in the tip holder and no tip in place | Contact the seller to repair the tip into place |
| Water tank leakage | Damaged seal ring | Contact the seller to repair or replace the seal ring |
| Water tank leakage | Low battery mode | Full charge |
| Failure to start | Low battery mode | Full charge |
| Full charge | Low battery mode | Cleaning of filter screen |

If the above solutions do not work, please contact the sales store or the designated after-sales service center.

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WARRANTY CONDITIONS:

The plant guarantees normal operation of the product within 12 months from the date of the sale, provided that the consumer complies with the rules of operation and care provided by this Manual. Service life of the product is 5 years. In order to avoid any misunderstanding we kindly ask you to carefully study the Owner's Manual, the terms of warranty, and check the correctness of filling the Warranty Card. The Warranty Card is valid only if the following information is correctly and clearly stated: model, serial number of the product, date of sale, clear seals of the company-seller, signature of the buyer. The model and serial number of the product must correspond to those specified in the Warranty Card. If these terms are violated, and if the information specified in the Warranty Card is changed, erased or rewritten, the Warranty Card is considered invalid. Setup and installation (assembly, connection, etc.) of the appliance are described in the attached documentation; they can be carried out both by the user, and by specialists having the necessary qualification sent by SC or selling companies (on a paid basis). In this case, the person (the institution) installing the appliance shall be responsible for the correctness and quality of installation (setup). Please pay attention to the importance of the correct installation of the appliance, both for its reliable operation, and for obtaining warranty and free service. Require the installation specialist to enter all the necessary information about the installation in the Warranty Card. In case of failure of the appliance during the warranty period by the fault of the manufacturer, the owner is entitled to free warranty repair upon presentation of the correctly filled Warranty Card together with the appliance in the factory configuration packed in the original packaging to the warranty workshop or the place of its purchase. Satisfaction of the customer's claim due to the fault of the manufacturer shall be made in accordance with the law of Ukraine «On Protection of Consumers' Rights». In the case of warranty repair the warranty period shall be extended for the period of repair and shipment.

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