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
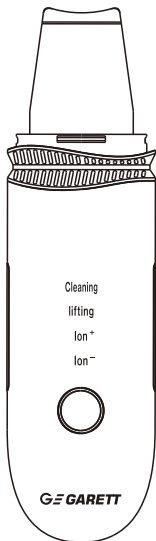
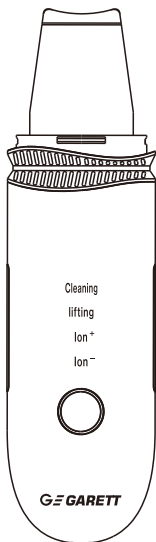
Beauty Sonic Scrub


Safety issues

1. Please read this manual carefully before using this product and keep it for future reference.
2. The product is charged with the charging kit provided.
3. It is forbidden to use this product while charging.
4. The product is not waterproof; do not submerge it in water.
5. Do not unplug the device from the outlet or use it with wet hands.
6. To avoid accidents, do not touch the metal part of the plug.
7. Do not destroy, modify, bend or twist the power cord - it may damage it.
8. Avoid getting the charging port and power cord dirty.
9. Charging should be stopped if the power plug or cord is damaged or becomes hot during charging.
10. Clean the device to avoid gathering dust and dirt.
11. Do not store the product in the bathroom, humid places, at high temperatures, as this may cause electric shock or fire.
12. Do not destroy the product, throw it or modify it. In the event of damage, stop using it.
13. Do not use the device around irritated, allergic or inflamed skin. In this case, please contact your doctor.
14. Children should not be allowed to use it or play with it.
15. Pregnant women or women during the menstrual cycle should not use the product without first consulting a doctor.

Product description

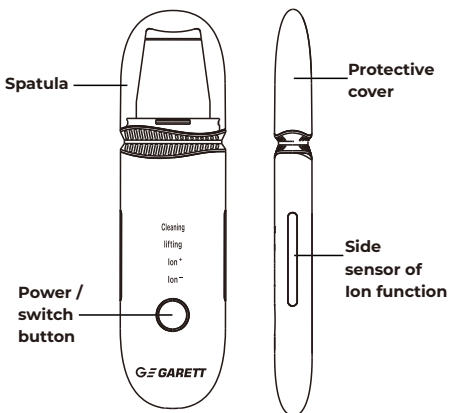Charging cable

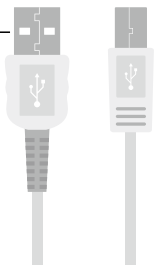Instructions for use

1. Make sure the device is charged before using it for the first time.
2. If there is no reaction after pressing the power button, connect the device to charging.
3. Connect the small end of the charging cable to the charging port and then to the USB adapter. While charging, the device shows the current battery level.
4. When the device is discharged, it takes about 2 hours to charge the battery fully.
5. To turn on the device, press and hold the power button for about 1.5 seconds - you will hear a sound notifying that the device is turned on.
6. Then select the operating mode by pressing the power / switch button. The change of the mode is indicated by a sound signal and the lighting of the current operating mode of the device.

Mode 1 – Cleansing the skin:

- Apply a cleanser or clean water onto the washed face and start the treatment.
- Move along the facial contours to be cleaned slowly.
- Remember to slope the spatula down against the pore direction.
- The treatment should take about 10 minutes.
When finished, cleanse your skin with clean water.


Mode 2 – Lifting:

- For basic use, clean your face with a steamer or apply a warm towel for 5 minutes.
- Apply a nourishing essence to your face and start the treatment.
- To remove blackheads, apply the balm for 10 minutes.
- Move along the facial contours to be cleaned slowly.
- Remember to slope the spatula down against the pore direction.
- The treatment should take about 10 minutes.
When finished, cleanse your skin with clean water.


Mode 3/4 – Ion Therapy (Positive Ion+ / Negative Ion-):

- For basic use, clean your face with a steamer or apply a warm towel for 5 minutes.
- Apply a nourishing essence to your face and start the treatment.
- Move along the facial contours to be cleaned slowly, starting from the bottom and moving upwards one by one.
- Remember to slope the spatula down against the pore direction.
- The treatment should take about 10 minutes. When finished, cleanse the skin, apply a cold mask, then keep your face in the steamer for 10 minutes or replace it with a warm towel. After the treatment, clean your face with clean water.

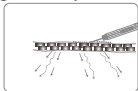Tips for user

Directions for the scapula movements during the procedure

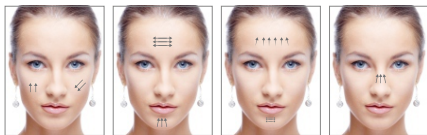1. When using the product, please keep your hand partially covering the side detectors for the correct operation of the 3/4 mode.
2. When using the product, the spatula angle should be about 45 degrees against the pore direction.
3. The product is a special cosmetic peeling device; do not use it for other purposes.
4. Do not use on areas with skin problems, after plastic surgery, sunburn, on irritated, allergic or inflamed skin.
5. It is recommended to use the device 2-3 times a week.
6. Before starting and after finishing the treatment, clean the spatula to extend the service life of the device.

Repair and maintenance

1. In the event of damage or failure, please stop using this product and contact the Garrett service. Any attempts to repair it on your own may damage the device and will invalidate the warranty.
2. Clean the device only with a soft, slightly damp cloth. Do not use chemicals or detergents for this purpose.

Technical problems

Symptoms	Possible reasons	Solutions
No vibration	Incorrect spatula position.	Check whether the spatula position has not changed, you can try to slightly correct it.
	Poor function or complete failure of the vibration module.	Contact the service.
It does not start	Uncharged battery.	Connect the device to charging and wait 2 hours.
	Battery damaged due to long period of disuse.	Contact the service.
	Damage caused by ingress of water into the interior of the device.	Contact the service.
Ion function does not work	No hand contact with side detectors.	During the procedure, hold the device in your hand so that it partially covers the side detectors.
	Defective button / switch.	Contact the service.
	Damage caused by ingress of water into the interior of the device.	Contact the service.

In case of further problems, please contact our Technical Support Department: serwis@garett.pl

Technical specification

Product name: Garrett Beauty Sonic Scrub

Dimensions: 195 × 96 × 38 mm

Weight: 86 g

Voltage: 3,7 V

Battery: 550 mAh

Spatula operation frequency: 28 kHz

Note: It is not recommended to be used by pregnant women or during the menstrual cycle, for people prone to allergic reactions and with skin problems. These people should not use the product without first consulting a physician.

Warranty / Complaint

The product is covered by a 24-month manufacturer's warranty. The warranty does not cover: defects and damage caused by the user's fault or as a result of improper use (not in accordance with the intended use), mechanical damage and resulting from the use of accessories other than those included in the set, consumables. When complaining about the product, please complete the warranty card attached to the device and send it together with the proof of purchase to the service address provided on the warranty card.

Contact details available on the website:

www.garett.eu

Producer:

PP S.C. Paweł Knez, Piotr Wilczyński

ul. Targowa 18/1413

25-520 Kielce

www.garett.eu

Correct Disposal of this Product
(Waste Electrical & Electronic Equipment)

The marking on the product or in related texts indicates that it should not be disposed of with other household wastes at the end of its working life. To avoid harmful effects on the environment and human health due to uncontrolled waste disposal, please separate the product from other types of waste and recycle responsibly to promote the reuse of material resources as permanent practice. Household users may obtain information on where and how to recycle this product for the environment in a safe way from the retailer where they purchased the product or from local authorities. Business users should contact their supplier and check the terms of the purchase contract. The product should not be disposed of with other commercial waste.

Made in P.R.C.