NÁVOD K POUŽITÍ
NÁVOD NA POUŽITIE
INSTRUKCJA OBSŁUGI
INSTRUCTION MANUAL
HASZNÁLATÍ ÚTMUTATÓ

PEP 712

Osobní váha s tělesnou analýzou
Osobná váha s telesnou analýzou
Waga łazienkowa z analizą ciała
Bathroom body analysis scale
Személyi mérleg test analízissal

Gallet

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OSTRZEŻENIE

Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

GENERAL SAFETY INSTRUCTIONS

- Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- **The appliance is intended for use in households only and for similar use! It is not intended for use in health-care facilities and for commercial use!**
- The appliance can use by people (including children) whose physical, sensual or mental inability or lack of experience and knowledge prevents safe use of the appliance if they are supervised.
- **Never use the appliance if it malfunctions, if it fell to the ground and got damaged or if it fell into water. In these cases bring the appliance to a specialised service to have the safety and functions of the device checked.**
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- **We recommend: Do not step onto the scale with wet feet or on wet surface of the scale, there is a risk of slipping.**
- **Never immerse the scale into water (even their parts)!**
- If the scale was stored in lower temperatures, let it adjust to the room temperature.
- Use the scale in a horizontal position only and at places where they cannot knock over, in sufficient distance from heat sources (e.g. **stove, fireplace, heat radiators**), sources of increased air humidity (e.g. **sauna, bathroom, swimming pool**) and equipment with strong electromagnetic field (e.g. **microwave oven, radio, mobile phone**).
- Protect the personal scale from dust, chemicals, excessive humidity, toiletries, liquid cosmetic preparations, intensive temperature changes and do not expose them to direct sunlight.
- Do not step on the edge of the scale! The scale could overturn and cause injury!
- Do not place the scale on uneven surface or on surface covered with carpet. Unstable or soft surface under the scale can have negative effect on precision of weighing.
- In handling the scale, make sure injury of persons or damage of the appliance is prevented.
- Do not overload the appliance with weight higher than the maximum capacity! You will thus prevent damaging the appliance
- Do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.
- Handle the scale with care. Do not throw them or jump on them. The scale is durable, but improper handling can destroy their electronic sensors.
- Do not disassemble the scale and do not take any parts out.
- Do not put any objects on the surface of the scale during storing, they could be damaged.
- If you do not intend to use the scale for a long time, remove batteries from them.
- If the battery is leaking, remove it immediately; it can damage the scale.
- Do not use worn and new batteries together and do not use rechargeable accumulators.
- Dispose of a discharged batteries in a suitable way.
- Never use the appliance for any other purpose than for that described in these instructions for use!

- Keep in mind that the data read from this scale are approximate and inaccurate in comparison to results of actual medicinal analyses. The exact determination of biometrical characteristics can be done only by a doctor via medical methods.
- Keep batteries and scale out of reach of children and incapacitated people.
A person, who would swallow the batteries, must seek medical help immediately.
- Do not use during pregnancy, with symptom of edema and during dialysis treatment.
- All text in other languages and pictures shown on the packaging, or the product itself, are translated and explained at the end of this language mutation.
- The producer is not responsible for damage caused by improper use of the appliance (e.g. breaking the board, etc.) and its guarantee does not cover the appliance in the case of failure to comply with the safety warnings above.

DESCRIPTION OF THE CONTROLS

P1. Unit of scale

P2. Display


P3. Glass surface

P4. Electrodes

P5. SET / ON button **SET**

P6. UP ▲ and DOWN ▼ buttons for setting

P7. KG / LB / ST button for setting

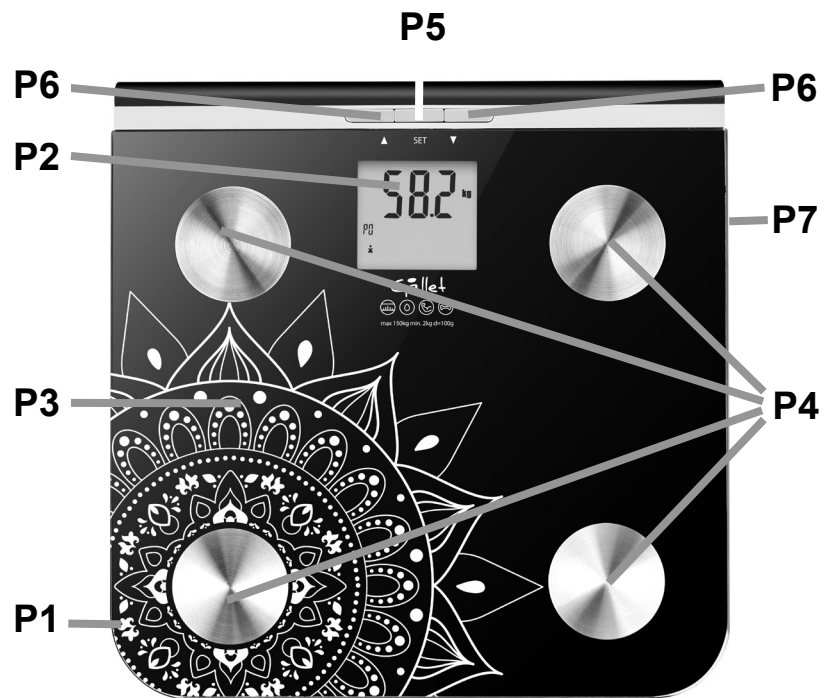LCD segmentum / icons illustrations:

P0 - P9 number of memory

MAM

WOMAN

FAT	measured value of fat (%)
TBW	measured value of water (%)
MUS	measured value of muscle (%)
BONE	bone weight (%)

SPECIFICATIONS

Measurement Range: 2.5 kg - 150 kg / st: 0.394 - 23.6 / 5.512 - 330.7 lb

Weight Measuring Division: 0.1 kg, 0.2 lb, 1/4 st

Water division: 0,1 %

Body fat division: 0,1 %

Age range for body fat: 10 - 85 years

Height range: 75 - 225 cm / 3ft 3.5in - 7 ft 2.5 in

Height unit: cm / ft-in

10 person memory

Recommended working environment: Temperature: 0 - 40 °C / 32 - 104 °F

Relative humidity: < 85 %

Choice of weighting modes: kilograms (kg) stones (st), pounds (lb)

BEFORE USE

Remove all the packing material and take the scale out. Remove all possible adhesion foils, stick-on labels or paper from the scale. Open the battery cover on the bottom of the scale, insert the battery with the correct polarity and close the cover in the opposite direction. Place the scale on flat, stable and dry surface. Turn on the scale by tapping your foot on the weighing surface and check the weight unit displayed on the display (**kg**). If you wish to set unit to **st** or **lb**, adjust the switch on the bottom of the scale **kg ↔ lb ↔ st (st:lb)** (or **cm ↔ FT;In**).

Then put the scale of firm and flat surface.

WEIGHT MEASURING

Why is knowledge of your weight important?

Overweight or underweight is an important indicator of your health. Overweight can increase the risk of some forms of cancer, heart diseases, diabetes, stroke, high blood pressure, arthritis, gall bladder disorders, apnoea and other diseases. Underweight brings with it risks, such as reduced immune system, osteoporosis, females can suffer from amenorrhoea which causes infertility etc. Obese people who want to lose weight and very thin people who want to gain weight should consult the way of treatment and diet with their doctor or dietician.

Why is it important to know what the value of fat and water in your organism is?

Obesity (excessive amount of body fat) is a known risk for people's health that is closely connected with serious diseases such as high blood-pressure, heart diseases, diabetes, osteoarthritis and some forms of cancer. There are other, less frequent symptoms, which may be caused by obesity, such as temporary suspension of breathing during sleep and tiredness during the day. An analysis of percentage ratio of body fat or your body composition is an effective way of assessing the level of your health. The contents of water in organism is a sign of good health. If the weight is not changing and the TBW % value (total contents of water in organism) is increasing, it means that the musculature in body is increasing (which is good for health). If the TBW % value is decreasing, it means decreasing of musculature and increasing of fat or it may indicate improper distribution of water in organism or some potential diseases relating to lack of water in the organism. In general, the more musculature, the higher the content of water in body.

For fat analysis it is essential:

- to wait approx. 15 minutes after getting up from bed because water in our body can be spread out.
- rest for a few hours (approx. 6–8) after an unusual exhausting activity.

The fat analysis will not provide useful information to children below 10 and adults over 85 years of age; persons with fever; after undergoing dialysis; persons with osteoporosis syndromes or persons who use cardiovascular medicaments; pregnant women; athletes who have a lot of musculature; because their results may be distorted.

How does the scale work?

The fat monitor uses the method of body composition analysis, BIA (short for bioelectric impedance analysis). BIA measures the resistance that is directly proportional to the fat in the body while the current flow is weak and harmless. Source electrodes let the current in on one side and detection electrodes detect this current on the other. For the purpose of analysis, data from the BIA appliance, such as sex, age and height, are subsequently inserted into a software program. The weight and body fat are calculated automatically. You will then find the percentage of fat and water in the body that is displayed in one of the ranges in the Reference table and find out which group you fall into. The accuracy of the scale depends on several variables of the user. We recommend making the weighing a routine, always weigh yourself approximately at the same time of the day and under same conditions.

Important notes:

- 1) This scale requires you to weigh yourself barefoot for the reason of correct conductivity during the body fat analysis.
- 2) Do not place the scale on an uneven surface or a surface covered with a carpet.
- 3) **Do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.**
- 4) Excessive drinking, eating, exercising, health conditions, used medicaments, women's menstruation cycle etc. influence the results of the analysis.

INSTRUCTIONS FOR USE

Your scale is fine electronic equipment. One of the things that are most difficult to weigh is human body as it is always in motion. If you want to obtain reliable results, always try to step with your feet at the same place on the scale and stand still (Fig. 1). If possible, weigh yourself at the same time of the day (best in the morning) and also on the same place to get comparable results. Remember that excessive clothes can increase your weight as well as eating and drinking before weighing. Short-term weight fluctuation is generally attributed to loss of liquids.

A) Regular weighing mode

You can use this scale only for weight measuring. This is basic mode without body fat measuring.

Turn the scale on by stepping on the weighing surface. After a moment, the measured weight flashes and the value remains displayed for several seconds. **Then the scale are switched off automatically.**

B) Weighing and total body analysis mode

Setting (change) of personal data for the body fat and water analysis

Before the first analysis, it is necessary to put in personal data (sex, height and age) which will be saved to the memory. It is possible to save data of 10 persons in the scale's memory.

The input data can be set and gradually changed in the following range:

Memory	Gender	Height	Age
0 – 9	male / female	75 – 225 cm	10 – 85

1. Press the **SET** button Number of memory will start the flash on the display (e.g. 1).
2. Select the number of the memory (with button ▲ or ▼) to which you want to save the data for the person to be weighed.
Confirm by pressing the **SET** button.
3. Set your sex with the ▲ or ▼ button:

♂ / female ♀

Confirm by pressing the **SET** button.

4. Set your age with the ▲ or ▼ button (Age)

Confirm by pressing the **SET** button.

5. Set your height with the ▲ or ▼ button (you can change **cm** ↔ **FT;In**).

Confirm by pressing the **SET** button.

Now all the personal data is saved (changed). Wait till 0.0 appears. Stand on the scale and wait for results.

The following data will start to show after a few seconds:

1. – Weight in set units
2. – Body fat contents in % (FAT) + analysis bar (- o + ++)
3. – Body water contents in % (TBW)
4. – Musculature contents in % (MUS)
5. – Bone weight contents in % (BONE)

The results will display in a cycle. After three repetitions, the scale will automatically turn off.

Notice

There must not be a contact between feet, calves or thighs. Otherwise the measurement could be inaccurate.

For a repeated analysis of body characteristics (after saving the input data of the person to the memory – e.g. 1), proceed as follows:

- 1) Press and hold the **SET** button (if you were using the scale in the regular weighing mode) or turn on the scale with the **SET** button.
- 2) Choose the number of memory with the ▲ or ▼ button where your setting is saved.
- 3) Do not press anything after selecting the respective number of memory (wait till 0.0 appears). Step on the scale and your weight and bodily characteristics will begin to be analyzed.

Age	Women		Man		
	Body Fat	Body Water	Body Fat	Body Water	
<= 30	4,0 - 20,5	66,0 - 54,7	4,0 - 15,5	66,0 - 58,1	Underweight
	20,6 - 25,0	54,6 - 51,6	15,6 - 20,0	58,0 - 55,0	Normal
	25,1 - 30,5	51,5 - 47,8	20,1 - 24,5	54,9 - 51,9	Overweight
	30,6 - 50,0	47,7 - 34,4	24,6 - 50,0	51,8 - 34,4	Obese
> = 30	4,0 - 25,0	66,0 - 51,6	4,0 - 19,5	66,0 - 55,4	Underweight
	25,1 - 30,0	51,5 - 48,1	19,6 - 24,0	55,3 - 52,3	Normal
	30,1 - 35,0	48,0 - 44,7	24,1 - 28,5	52,2 - 49,2	Overweight
	35,1 - 50,0	44,6 - 34,4	28,6 - 50,0	49,1 - 34,4	Obese

Underweight indicator (on scale): -
Normal weight indicator (on scale): o
Overweight indicator (on scale): +
Obesity indicator (on scale): ++

For the reason of differences between people, this table is only informative.
The results of measuring bones' weight and muscular proportions are simply tentative.
 Healthy man has approximately 40% of muscular proportion and healthy woman has 30% according to an average ratio. Weight of bones depends on many relevant factors. Always consult your attending physician if any measuring results are abnormal.

Note

The scale is switched on automatically if you do not step on it within approx. 20 seconds after „0.0“ is displayed.

ERROR MESSAGES AND THEIR SYMBOLS

Symbols displayed	Interpretation
Lo	Battery is dead
Err	Very low / high % of fat
EEEE	The measured weight is higher than maximum allowed value of the scale.

TROUBLESHOOTING OF POSSIBLE PROBLEMS

If you have problems with weight:

- Check whether the batteries are inserted properly.
- Check whether you selected the proper weight unit.
- Check whether the scale is on a flat floor and they do not touch the wall.
- If nothing appears on the scale when you want to switch them on by tapping on them or if „Lo“ appears only, replace the battery.
- If „EEEE“ appears on the display, the scale was overloaded.
- If unusual signs are displayed, remove the insert the battery again.
- If signs appear, such as if the display is flashing or if there is a faulty sign, try to move the scale away from the influence of a source of interference or ensure switching off the source for the time of using the scale.

STORING

After using the scale, put it into the horizontal position (not vertically), out of reach of children and incapacitated people.

MAINTENANCE


Do not use coarse and aggressive cleaning agents (e.g. **sharp objects, scrapers, diluting agents or other solvents**). Clean the board by wiping with a soft wet cloth. Make sure that water does not get inside the appliance.

TECHNICAL SPECIFICATIONS

- Electronic glass scale with body fat, water, bone weight and muscle percentage
- Elegant inverse LCD display
- Auto on – immediately show you weight;
- 6 mm tempered glass surface
- Using 5 parameters for measuring of body fat – height, weight, age, sex and biological resistant
- Auto off after 20 sec.; auto zero; Overload indicator; Low battery indicator
- Batteries operated 2 x AAA (included)
- Size: 30.5 x 30.5 x 2 cm; Weight: 2 kg
- Protection class III, CE Certification

The manufacturer reserves the right to change the technical specifications and accessories for the respective models. Product image are for illustrative purposes only.

Old electrical appliances, used batteries and accumulators disposal

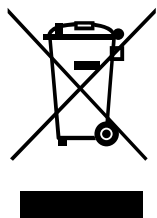This symbol appearing on the product, on the product accessories or on the product packing means that the product must not be disposed as household waste. When the product/ battery durability is over, please, deliver the product or battery (if it is enclosed) to the respective collection point, where the electrical appliances or batteries will be recycled. The places, where the used electrical appliances are collected, exist in the European Union and in other European countries as well. By proper disposal of the product you can prevent possible negative impact on environment and human health, which might otherwise occur as a consequence of improper manipulation with the product or battery/ accumulator. Recycling of materials contributes to protection of natural resources. Therefore, please, do not throw the old electrical appliances and batteries/ accumulators in the household waste. Information, where it is possible to leave the old electrical appliances for free, is provided at your local authority, at the store where you have bought the product. Information, where you can leave the batteries and accumulators for free, is provided to you at the store, at your local authority.

HOUSEHOLD USE ONLY.

DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS.

AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements.

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