DK

NO

SE

FI

GB

DE

PL

ET

ES

IT

NL

FR

Brugsanvisning · Bruksanvisning · Bruksanvisning  
 Käyttöohje · Instruction manual · Gebrauchsanweisung  
 Podręcznik użytkownika · Kasutusjuhend  
 Manual de instrucciones · Manuale di istruzioni  
 Gebruiksaanwijzing · Manuel d'instructions

# Ledningsfri massagepistol

DK

## INTRODUKTION

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# Cordless Massage Gun

## INTRODUCTION

To get the most out of your new product, please read this manual before use. Please also save the manual in case you need to refer to it at a later date.

## TECHNICAL DATA

Battery:	Lithium battery, 1800mAh
Charging cable:	Type-C, DC 5 V 1 A
Charging time:	3.5 hours
Rated usage time:	15 minutes
Operating time after full charge:	Up to 1.5 hours
Speed:	6 settings, 2400-3200 bpm
Includes	6 massage heads

Do not use this appliance in a bathtub, shower, or water-filled bowl.

Class III product: Designed to be supplied from a separated extra low voltage supply (SELV). The voltage from a SELV supply is low enough that a person can normally come into contact with it without the risk of electric shock.

## SPECIAL SAFETY INSTRUCTIONS

This appliance must not be used by persons, including children, with a physical, sensory or mental disability or anyone without experience or knowledge unless they are supervised or have been instructed how to use the appliance.

Children should be supervised to ensure that they do not play with the appliance.

If you are pregnant, have a pacemaker or have diabetes or any other health issues, you should ask your doctor for advice before using the massage gun.

If you have persistent muscle or joint pain that does not disappear after using the massage gun, you should ask your doctor for advice.

Do not use the massage gun on areas of skin with open wounds, or if the skin is discoloured, burned, inflamed or affected by eczema.

Do not use the massage gun on wet or damp skin. If necessary, use the massage gun over your clothes.

Do not use the massage gun on the head or bones, e.g. on the knee or elbow.

Keep hair and fingers away from the massage head to avoid them becoming trapped.

The massage gun must not be used for more than 15 minutes at a time.

Check regularly that the massage gun and charger are in good working order. If this is not the case, the massage gun must not be used.


Do not attempt to repair the massage gun yourself.

The massage gun and the charger must not be exposed to heat, direct sunlight, moisture, sharp edges and the like.

Always turn the massage gun off when it is not in use or when cleaning it.

The massage gun is solely intended for domestic use and is not suitable for commercial use. Always use the massage gun in accordance with the instructions.

## MAIN COMPONENTS

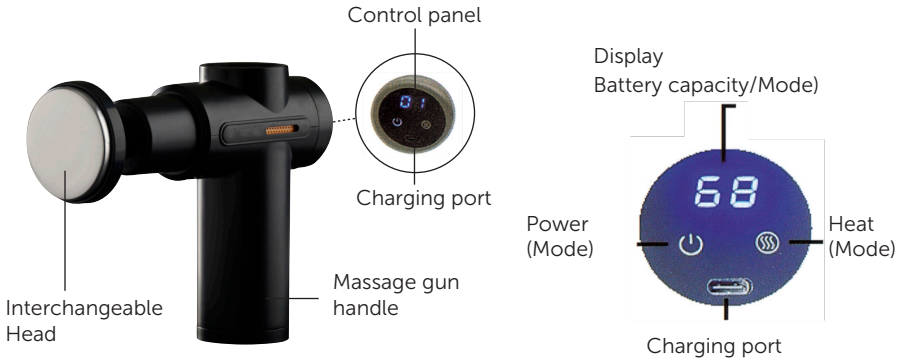Selected Mode Display

Heat Mode OFF

Heat Mode ON

Hot Head  
(To soothes stiff, aching muscles)

Cold Head  
(For joint swelling and muscle inflammation)

Round Head  
(For large muscle groups)

Flat Head  
(For any parts of body)

Fork Head  
(For neck & spine)

Bullet Head  
(For joints)

## USE

Find yourself a suitable seated position. It is not recommended to use the massage gun when standing or when prone.

Choose a suitable massage head and plug it into the massage gun by gently pressing to secure it. To remove a head, gently pull to remove. Do not add or remove massaging heads whilst the massage gun is ON.

Speed setting	Beats per minute
1	2400
2	2600
3	2800
4	2900
5	3000
6	3200

To turn the massage gun ON/OFF press the Power (Mode) button.

Start at the lowest power setting and gently apply the massage gun to your body to start your massage.

Short press the Power (Mode) button to cycle through the available 6 power levels (where 1 is the gentlest and 6 is the most vigorous). Choose a power level where you feel comfortable. Excessive pressure can cause injury and may damage the device.

### Note:

The massage gun will automatically turn off after 15 minutes. This is the recommended massage time. Exceeding the recommended massage time may cause injury and/or damage to the massage gun.

If the massage gun is on and you press the massage head too hard, the massage gun will stop working and "Er" will show on the display. Please do not worry, this is a built-in feature to avoid damaging the motor.

## HOW TO USE HOT-HEAD ATTACHMENT

Ensure that the massage gun is turned OFF before attaching the Heating Head to the massage gun.

Turn the massage gun ON by long pressing the Power (Mode) button. Long press the Heat (Mode) button for 2 seconds. The Heating mode symbol will illuminate.

Place the massage gun on a suitable flat, dry, heat resistant surface and allow 3-4 minutes for the Heating Head to reach the appropriate temperature.

Press the Power (Mode) button to set the appropriate power level.

Apply the massage gun to your body to start your massage.

**Note:** Exercise caution to avoid body parts which are sensitive to heat to prevent burns.

When finished, hold the Power (Mode) button to turn off the massage gun. Place onto a suitable flat, dry heat resistant surface and allow to cool fully before adding or removing any accessories or recharging.

## HOW TO USE COLD HEAD ATTACHMENT

Before using the Cold Head attachment, detach it from the massage gun if necessary and place the head into a refrigerator for around 2 hours. Do not cool below 4°C.

Attach the Cold Head to the massage gun.

Turn the massage gun ON by long pressing the Power (Mode) button.

Start on the lowest power setting and slowly cycle through the power levels to find your comfort level.

Apply the massage gun to your body to start your massage.

**Note:** Exercise caution to avoid body parts which are sensitive to cold.

When finished, hold down the Power (Mode) button to turn the massage gun off.

## CHARGING THE MESSAGE GUN

**Before using the message gun for the first time, it needs to be charged fully (about 3-4 hours).**

Connect the USB Charging Cable to the message gun and a power source.

The battery capacity indicator will turn on when the battery is being charged and display the current battery level.

The battery capacity indicator will display '99', when the battery is fully charged.

The battery can be recharged at any time and at any level.

When the product is OFF, short press the Power (Mode) button to show the remaining battery capacity.

### Note:

Motor noises during operation is normal.

If the rated working time is over, the product will automatically turn off.

If you use this product continuously for over 30 minutes, the temperature control will automatically turn off the device. If this occurs, please allow for a 30 minute break before resuming.

## CLEANING AND MAINTENANCE

Before cleaning, please make sure the message gun has been turned off and the adapter has been disconnected from the product.

Use slightly damp cloth to clean the surface of product and control panel.

### Environmental information

Electric and electronic equipment (EEE) contains materials that could be dangerous and harmful to human health and the environment if the electrical and electronic waste (WEEE) is not properly disposed of. Products (and batteries if included) marked with the crossed-out waste bin are electrical and electronic equipment. The crossed-out waste bin symbolises that the device must not be disposed of at the end of its service life together with household waste but must be disposed of separately.

## EU DECLARATION OF CONFORMITY

<https://manuals.media/details/?id=100009>

---

## SERVICE CENTRE

Note: Please provide the product model number in connection with all inquiries. The model number is shown on the front of this manual and on the product rating plate. [www.schou.com](http://www.schou.com)

---

Manufactured in China  
Manufacturer:  
Schou Company A/S  
Nordager 31  
DK-6000 Kolding

All rights reserved. The content of this manual may not be reproduced, either in full or in part, in any way by electronic or mechanical means, e.g. photocopying or publication, translated or saved in an information storage and retrieval system without written permission from Schou Company A/S.

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